



Protein Phosphatase Cdc25C Fluorometric Assay Kit

User's Manual

For Research Use Only, Not for use in diagnostic procedures

Fluorometric Assay Kit for Measuring Cdc25C Phosphatase Activity

CycLex Protein Phosphatase Cdc25C Fluorometric Assay Kit

100 Assays

Cat# CY-1354

Intended Use..... 1

Storage..... 1

Introduction 2

Principle of the Assay..... 2

Materials Provided 3

Materials Required but not Provided..... 3

Precautions and Recommendations..... 3-4

Detailed Protocol..... 4-6

Evaluation of Results 7-9

Troubleshooting 10

Reagent Stability 10

References..... 11

Related Products..... 11-12

Intended Use

The CycLex Research product **Protein Phosphatase Cdc25C Fluorometric Assay Kit** Users Manual is a fluorometric and non-radioactive assay designed to measure the activity of Cdc25C protein phosphatase. This 96-well assay is useful for screening inhibitors and modulators of Cdc25C activity in HTS. The kit includes all necessary components, including recombinant, human Cdc25C (catalytic domain), for use in preinvestigational drug discovery assays.

This assay kit is for research use only and not for use in human, diagnostic, or therapeutic procedures.

Storage

- Upon receipt, store the kit at -70°C.
- Don't expose reagents to excessive light.
- **AVOID REPEATED FREEZE THAW CYCLES OF "Recombinant Cdc25C"!**



Introduction

Activation of cyclin-dependent kinases in higher eukaryotic cells can be achieved through dephosphorylation of two conserved residues, Thr14 and Tyr15 (1) by members of the Cdc25 phosphatase family, Cdc25A, Cdc25B and Cdc25C. The Cdc25 dual-specificity phosphatases control progression through the eukaryotic cell division cycle by activating cyclin-dependent kinases. Cdc25A plays an important role at the G1/S-phase transition. Cdc25A degradation during mitotic exit and in early G1 is mediated by the anaphase-promoting complex or cyclosome (APC/C)(Cdh1) ligase, and that a KEN-box motif in the N-terminus of the protein is required for its targeted degradation (2, 3). Cdc25B undergoes activation during S-phase and plays a role in activating the mitotic kinase Cdk1/cyclin B in the cytoplasm (4-6). Active Cdk1/cyclin B then phosphorylates and activates Cdc25C leading to a positive feedback mechanism and to entry into mitosis. In addition to their essential function in the normal cell cycle, cdc25 phosphatases are involved in the checkpoint-induced control of cell cycle progression (7-9). Finally, Cdc25A and B phosphatases have been shown to possess an important oncogenic potential and to be overexpressed in a variety of cancers and cancer cell lines (4, 10).

Principle of the Assay

The **Protein Phosphatase Cdc25C Fluorometric Assay Kit** is based on an exclusive fluorescence substrate, OMFP (3-o-methylfluorescein phosphate). This homogenous assay kit is sensitive and convenient. This method of measurement should raise the efficiency of inhibitor screening and biochemical analysis of this enzyme.

Summary of Procedure

Mix 40 μ L of Assay mixture and 5 μ L of test compound in the wells



Add 5 μ L of Recombinant Cdc25C



Incubate for 15 min at room temp.

Add 25 μ L of Stop Solution



Measure fluorescence at 510-540 nm emission / 482-502 nm excitation



Protein Phosphatase Cdc25C Fluorometric Assay Kit

User's Manual

For Research Use Only, Not for use in diagnostic procedures

Materials Provided

All samples and standards should be assayed in duplicate. The following components are supplied and are sufficient for one hundred assays.

Components of Kit

| Components | Quantity | Storage |
|--|-------------------|-------------|
| ① 10X Assay Buffer | 600 μ L x 1 | Below -20°C |
| ② 10X 3-o-methyl fluorescein phosphate (OMFP) | 550 μ L x 1 | Below -20°C |
| ③ Recombinant Cdc25C, Human | 550 μ L x 1 | -70°C |
| ④ 100X Cdc25 Inhibitor : 100 μ M Na ₃ VO ₄ in DW | 100 μ L x 1 | Below -20°C |
| ⑤ Stop Solution | 1,300 μ L x 2 | Below -20°C |
| ⑥ Instruction Manual | 1 | Room temp. |

Materials Required but not Provided

- **Microtiter plate suitable for use with a fluorometric plate reader**
- **Fluorometric plate reader or microtiter plate fluorometer:** Use a fluorescence microplate reader equipped with appropriate filters. OMFP has excitation/emission maxima of approximately 485/525 nm. We have found that standard filters for blue-fluorescent dyes (e.g., excitation = 485 \pm 12.5 nm, emission = 525 \pm 20 nm) can be used to detect OMFP.
- **Pipettors:** 2-20 μ L, 20-200 μ L and 200-1000 μ L precision pipettors with disposable tips
- **Multi-channel pipette**
- **Microtiter plate shaker**
- **Distilled water (DW)** or equivalent high quality water
- **Microcentrifuge and tubes** for sample preparation
- **Reagent reservoirs**
- **Ice bucket** to keep reagents cold until use

Precautions and Recommendations

- Upon receipt, store the kit at -70°C.
- Do not expose reagents to excessive light.
- Do not use kit components beyond the indicated kit expiration date.
- Rinse all detergent residue from glassware.
- Use deionized water of the highest quality.
- Do not mix reagents from different kits.
- Do not mouth pipette or ingest any of the reagents.
- Do not smoke, eat, or drink when performing the assay or in areas where samples or reagents are handled.



Protein Phosphatase Cdc25C Fluorometric Assay Kit

User's Manual

For Research Use Only, Not for use in diagnostic procedures

NOTE: THE FOLLOWING PROCEDURES ARE INTENDED ONLY AS A GUIDELINE. THE OPTIMAL EXPERIMENTAL CONDITIONS WILL VARY DEPENDING ON THE PARAMETERS BEING INVESTIGATED, AND MUST BE DETERMINED BY THE INDIVIDUAL USER.

Detailed Protocol

Preparation of Reagents

Thaw the reagents at room temperature except “Recombinant Cdc25C” and keep all reagents on ice until use. Use them only after they are completely thawed and mixed.

1. Prepare **10X Cdc25 Inhibitor** by adding 5 μL of the ④ **100X Cdc25 Inhibitor** (provided) to 45 μL of distilled (deionized) water. Mix well.

Discard any unused 10X Cdc25 Inhibitor after use.

2. Prepare **Assay Mixture** by adding 5 μL of the ①10X Assay Buffer (provided) and 5 μL of the 10X Fluoro-Phospho-Substrate (provided) to 30 μL of distilled (deionized) water per one assay. Mix well.

Assay Mixture

| Assay reagents | 1 assay | 8 assays | 16 assays | 32 assays | 48 assays |
|--------------------------------------|------------------------------------|-------------------------------------|-------------------------------------|---------------------------------------|---------------------------------------|
| Distilled water | 30 μL | 240 μL | 480 μL | 960 μL | 1,440 μL |
| ①10X Assay Buffer | 5 μL | 40 μL | 80 μL | 160 μL | 240 μL |
| ②10X OMFP | 5 μL | 40 μL | 80 μL | 160 μL | 240 μL |
| Total volume of Assay Mixture | 40 μL | 320 μL | 640 μL | 1,240 μL | 1,920 μL |



Assay Procedure

In order to estimate the inhibitory effect on Cdc25C activity by the test compounds correctly, it is necessary to conduct the control experiment of “**Vehicle control**” at least once for every experiment and “**Inhibitor control**” at least once for the first experiment, in addition to “**Test sample**” as indicated in the Table.1 (below). When test chemicals cause an inhibitory effect on Cdc25C activity, the level of increase of fluorescence intensity is weakened as compared with “**Vehicle control**”. The increase in fluorescence intensity is not observed in “**Inhibitor control**”.

1. Following Table.1 below, first, add “**Assay mixture**” to microtiter plate wells. Second, add “**Test Compound**” or “**Vehicle of Test Compounds**” or “**10X Cdc25 Inhibitor**” to each well of the microtiter plate and mix well.

Table.1: Reaction mixture

| Assay reagents | Test Sample | Vehicle Control | Inhibitor Control | No Enzyme Control |
|---|-----------------------------|-----------------------------|-----------------------------|-----------------------------|
| Assay Mixture | 40 μ L | 40 μ L | 40 μ L | 40 μ L |
| Test Compound | 5 μ L | - | - | - |
| Vehicle of Test Compounds | - | 5 μ L | - | 5 μ L |
| 10X Cdc25 Inhibitor* | - | - | 5 μ L | - |
| ③Recombinant Cdc25C (0.2 μ g/ μ L) | 5 μ L | 5 μ L | 5 μ L | - |
| Distilled water | - | - | - | 5 μ L |
| Total Volume of the Reaction mixture | 50 μL | 50 μL | 50 μL | 50 μL |

*10X Cdc25 Inhibitor (10 μ M Na_3VO_4): See Page 4, section Preparation of Reagents

2. Initiate reactions by adding 5 μ L of “**③Recombinant Cdc25C**” or distilled water to each well and mixing thoroughly at room temperature.
3. Incubate for 15 min or desired length of time at room temperature.
4. Add 25 μ L of “**⑤Stop Solution**” to each well of the microtiter plate, and mix thoroughly.
5. Measure fluorescence intensity using a microtiter plate fluorometer with excitation at 482-502 nm and emission at 510-540 nm.
6. The efficacy of the Test compound is the difference in fluorescence intensity between “**Vehicle control**” and “**Test sample**”.

Note: If necessary, it is possible to store the microtiter plate after adding “**④Stop Solution**” for a few hours at 4°C. The microtiter plate must be sealed to prevent evaporation and kept from excessive light.



Protein Phosphatase Cdc25C Fluorometric Assay Kit

User's Manual

For Research Use Only, Not for use in diagnostic procedures

Alternate procedure

- 1'. Following Table.1 above, first, add **"Assay mixture"** to microtiter plate wells. Second, add **"Test Compound"** or **"Vehicle of Test Compounds"** or **"10X Cdc25 Inhibitor"** to each well of the microtiter plate and mix well.
- 2'. Initiate reactions by adding 5 μ L of **"Recombinant Cdc25C"** or distilled water to each well and mixing thoroughly at room temperature.
- 3'. Read fluorescence intensity for 20 to 30 minutes at 1 to 2 minute intervals using microtiter plate fluorometer with excitation at 482-502 nm and emission at 510-540 nm.
- 4'. Measure and calculate the rate of reaction while the reaction velocity remains constant.

Caution and Significance

- All samples and "Recombinant Cdc25C" should be assayed in duplicate.
- Use of a microtiter plate shaker is recommended for complete mixing.
- If the test compounds or samples themselves emit fluorescence at excitation wavelength: 482-502 nm and fluorescence wavelength: 510-540 nm, the test assay cannot be evaluated correctly.



Evaluation of Results

Analysis of Inhibitor Effect #1


% Intensity

1. Run reactions with test compounds and Vehicle as described in the **Detailed Protocol**.
2. Subtract fluorescence intensity of "No Enzyme Control" from all experimental samples (Test Samples and Vehicle Control).
3. Calculate the % Intensity:

$$\% \text{ Intensity} = \frac{\text{Fluorescence Intensity of Test Sample}}{\text{Fluorescence Intensity of Vehicle Control}} \times 100$$

Note: This % Intensity is a rough value of enzyme activity or inhibition. For greater accuracy, plot a standard curve of Cdc25C for each new set of reactions and estimate the % Activity (see below).

Fig.1 Cdc25C Inhibition Curve by SOV(Na_3VO_4 ,sodium orthovanadate)

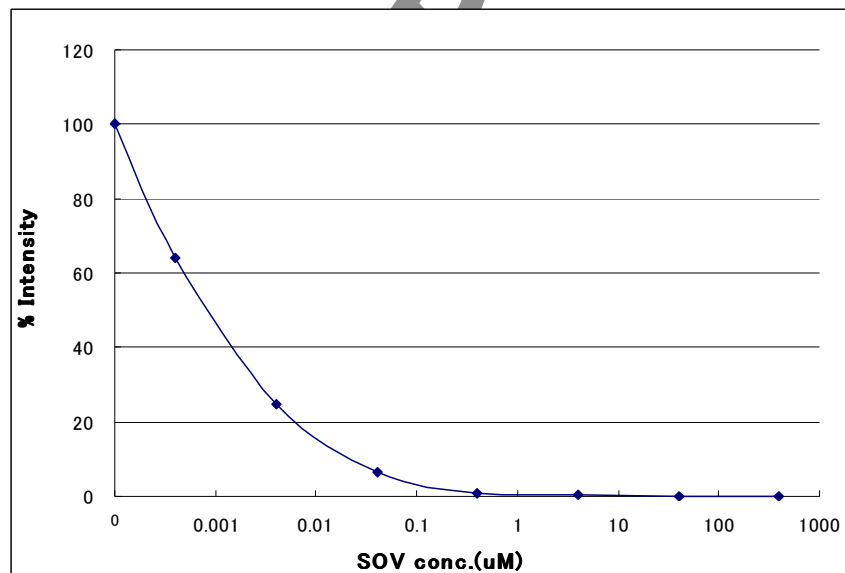




Analysis of Enzyme Activity

Cdc25C Standard Curve and % Activity

1. Dilute the ①10X Assay buffer 1:10 with distilled water to make 1X Assay Buffer.
2. Make serial dilutions of Recombinant Cdc25C with 1X Assay Buffer (ex. 100%, 50%, 25%, 12.5%, 6.25%, 3.13% and 0%).
3. Run reactions with Vehicle and serial dilutions of Recombinant Cdc25C as described in the Detailed Protocol.

1. Plot standard curve data (dose dependent curve data) as fluorescence intensity at 510-540 nm versus dose of Cdc25C (ng/assay)
2. Obtain a line-fit to the data using appropriate calculations.
3. Use the slope and Y-intercept to calculate the amount of Cdc25C activity for the experimental data.

Fig.2 Dose Dependency of Recombinant Cdc25C catalytic domain

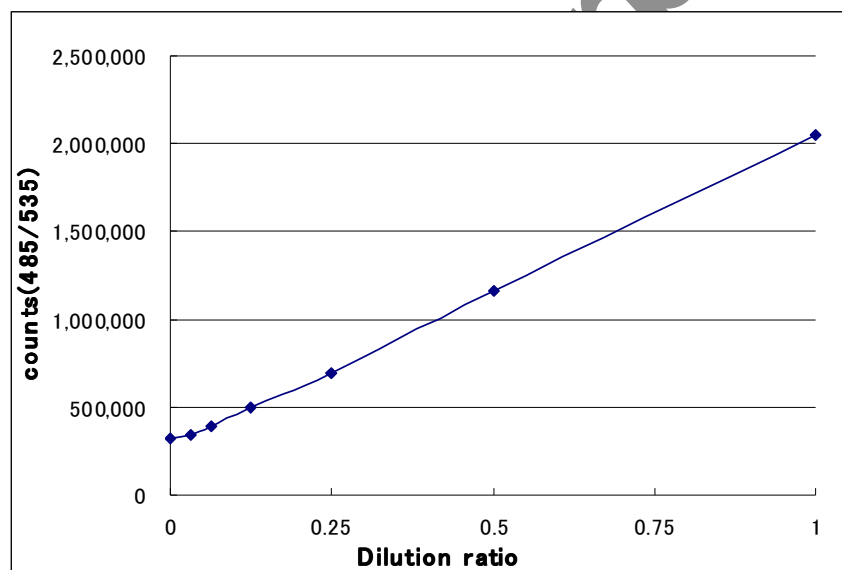



Analysis of Kinetics


Time Course Curve

1. Run reactions as described in the *Detailed Protocol*.
2. Subtract fluorescence intensity at the 0 time from all reaction time points.
3. Plot fluorescence intensity at 510-540 nm versus reaction time.
4. Determine the reaction time range in which the increase in fluorescence intensity at 510-540 nm is linear.
5. Calculate activity:

$$\text{Activity (reaction velocity)} = \frac{\text{Fluorescence Intensity of Test Sample}}{\text{Reaction time (min.)}}$$

Note: Usually, the linear range is from 0 to 30 min. This value is variable depending on reaction conditions and storage/handling of the Recombinant Cdc25C. Decreasing the amount of Recombinant Cdc25C in the assay may help to lengthen the time range.

Fig.3 Time Course Curve of Recombinant Cdc25C

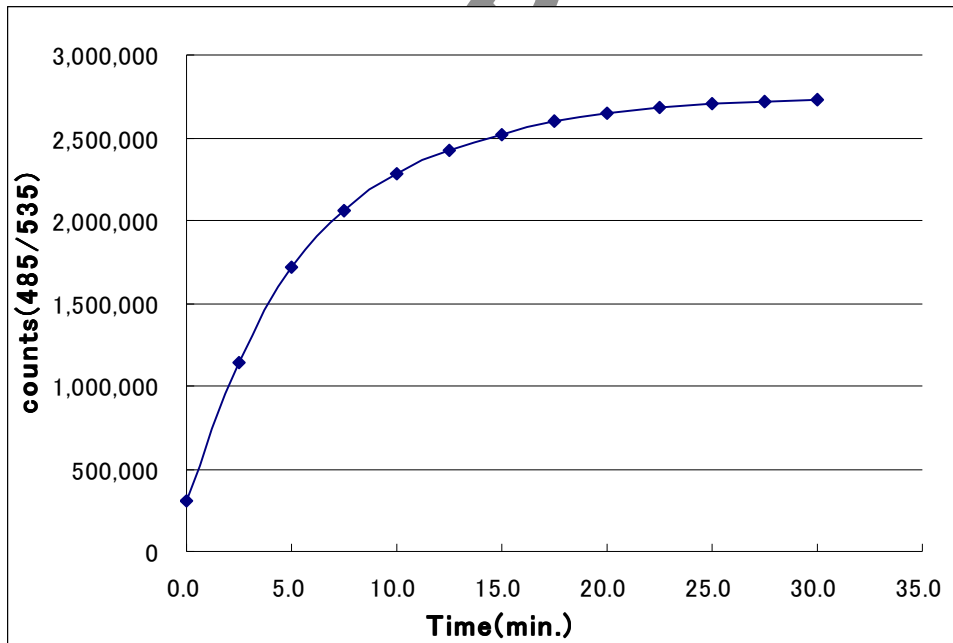



Troubleshooting

1. The Recombinant Cdc25C should be run in duplicate using the protocol described in the *Detailed Protocol*. Incubation times or temperatures significantly different from those specified may give erroneous results.
2. The reaction curve is nearly a straight line if the kinetics of the assay is of the first order. Variations in the protocol can lead to non-linearity of the curve, as can assay kinetics of other than first order. For a non-linear curve, point to point or quadratic curve fit methods should be used.
3. Poor duplicates, accompanied by elevated values for wells containing no sample, indicate inaccurate dispensing of assay reagents. If all instructions in the *Detailed Protocol* were followed accurately, such results indicate a need for multi-channel pipette maintenance.

Reagent Stability

All of the reagents included in the **Protein Phosphatase Cdc25C Fluorometric Assay Kit** have been tested for stability. Reagents should not be used beyond the stated expiration date. Upon receipt, all kit reagents, except Recombinant Cdc25C, should be stored at -20°C . Recombinant Cdc25C should be stored at -70°C . After use, return kit reagents to -20°C as soon as possible.

For research use only, not for use in human, diagnostic or therapeutic procedures



References

1. Morgan, D. O. (1997) *Annu. Rev. Cell Dev. Biol.* **13**, 261-291
2. Jinno, S., Suto, K., Nagata, A., Igarashi, M., Kanaoka, Y., Nojima, H., and Okayama, H. (1994) *EMBO J.* **13**, 1549-1556
3. Hoffmann, I., Draetta, G., and Karsenti, E. (1994) *EMBO J.* **13**, 4302-4310
4. Nagata, A., Igarashi, M., Jinno, S., Suto, K., and Okayama, H. (1991) *New Biol.* **3**, 959-968
5. Karlsson, C., Katich, S., Hagting, A., Hoffmann, I., and Pines, J. (1999) *J. Cell Biol.* **146**, 573-583
6. Hoffmann, I., Clarke, P. R., Marcote, M. J., Karsenti, E., and Draetta, G. (1993) *EMBO J.* **12**, 53-63
7. Furnari, B., Rhind, N., and Russell, P. (1997) *Science* **277**, 1495-1497
8. Sanchez, Y., Wong, C., Thoma, R. S., Richman, R., Wu, Z., Piwnica-Worms, H., and Elledge, S. J. (1997) *Science* **277**, 1497-1501
9. Peng, C.-Y., Graves, P. R., Thoma, R. S., Wu, Z., Shaw, A. S., and Piwnica-Worms, H. (1997) *Science* **277**, 1501-1505
10. Galaktionov, K., Lee, A. K., Eckstein, J., Draetta, G., Meckler, J., Loda, M., and Beach, D. (1995) *Science* **269**, 1575-1577

Related Products

- * CycLex Protein Tyrosine Phosphatase 1B (PTP1B) Fluorometric Assay Kit: Cat# CY-1350
- * CycLex T Cell Protein Tyrosine Phosphatase (TC-PTP) Fluorometric Assay Kit: Cat# CY-1351
- * CycLex Protein Phosphatase Cdc25A Fluorometric Assay Kit: Cat# CY-1352
- * CycLex Protein Phosphatase Cdc25B Fluorometric Assay Kit: Cat# CY-1353
- * CycLex Protein Phosphatase Cdc25C Fluorometric Assay Kit: Cat# CY-1354
- * CycLex Protein Phosphatase Cdc25 Combo Fluorometric Assay Kit: Cat# CY-1355
- * CycLex Protein Phosphatase Cdi1/KAP Fluorometric Assay Kit: Cat# CY-1356
- * CycLex Protein Phosphatase LMW-PTP/ACP1 Fluorometric Assay Kit: Cat# CY-1358
- * CycLex Protein Phosphatase DUSP1/MKP-1 Fluorometric Assay Kit: Cat# CY-1373

- * Protein Tyrosine Phosphatase PTPRA 1st Catalytic Domain: Cat# CY-E1301
- * Protein Tyrosine Phosphatase PTPRA 2nd Catalytic Domain: Cat# CY-E1302
- * Protein Tyrosine Phosphatase PTPRD 2nd Catalytic Domain: Cat# CY-E1307
- * Protein Tyrosine Phosphatase PTPRE 1st Catalytic Domain: Cat# CY-E1308
- * Protein Tyrosine Phosphatase PTPRF 1st Catalytic Domain: Cat# CY-E1310
- * Protein Tyrosine Phosphatase PTPRK 1st Catalytic Domain: Cat# CY-E1316
- * Protein Tyrosine Phosphatase PTPRQ: Cat# CY-E1323
- * Protein Tyrosine Phosphatase PTP4A2: Cat# CY-E1341
- * Recombinant Cdc25A (Catalytic domain): Cat# CY-E1352
- * Recombinant Cdc25B (Catalytic domain): Cat# CY-E1353
- * Recombinant Cdc25C (Catalytic domain): Cat# CY-E1354
- * Recombinant Cdi1/KAP: Cat# CY-E1356
- * Protein Phosphatase PP5: Cat# CY-E1359



Protein Phosphatase Cdc25C Fluorometric Assay Kit

User's Manual

For Research Use Only, Not for use in diagnostic procedures

- * Protein Tyrosine Phosphatase PTPN3/PTPH1: Cat# CY-E1360
- * Protein Tyrosine Phosphatase PTPN6/SHP-1: Cat# CY-E1363
- * Protein Tyrosine Phosphatase PTPN7/HePTP: Cat# CY-E1364
- * Protein Tyrosine Phosphatase PTPN8/PTPN22: Cat# CY-E1365
- * Protein Tyrosine Phosphatase PTPN9/MEG2: Cat# CY-E1366
- * Protein Tyrosine Phosphatase PTPN11/SHP-2: Cat# CY-E1367
- * Protein Tyrosine Phosphatase PTPN12/PTP-PEST: Cat# CY-E1368
- * Protein Tyrosine Phosphatase PTPN13/FAP-1: Cat# CY-E1369
- * Protein Tyrosine Phosphatase PTPN14/PEZ: Cat# CY-E1370
- * Protein Tyrosine Phosphatase PTPN21/PTPD1: Cat# CY-E1372
- * Protein Phosphatase DUSP1/MKP-1: Cat# CY-E1373

PRODUCED BY

CycLex Co., Ltd.
1063-103 Terasawaoka
Ina, Nagano 396-0002
Japan
Fax: +81-265-76-7618
e-mail: info@cyclex.co.jp
URL: <http://www.cyclex.co.jp>

CycLex/CircuLex products are supplied for research use only. CycLex/CircuLex products and components thereof may not be resold, modified for resale, or used to manufacture commercial products without prior written approval from CycLex Co., Ltd.. To inquire about licensing for such commercial use, please contact us via email.